

Adapter and Trigger Sprayer

Catalog No. F11620-0050

The SP Bel-Art Adapter has one end to fit the opening on all 500ml and 1000ml SP Bel-Art wash bottles and another end to fit the SP Bel-Art Trigger Sprayer. The adapter and trigger sprayer are compatible with a wide range of chemicals. However, because the sprayer can be adjusted to produce a fine mist, it is not recommended for use near open flames or other sources of combustion. Directions for affixing the adapter and trigger sprayer are provided below.

Directions for Assembly:

1. To attach adapter to a 1000ml SP Bel-Art wash bottle, skip to Step #7
2. To attach adapter to a 500ml SP Bel-Art wash bottle, you will need to cut the tubing as described below
3. Grasp the tubing and pull it away from the trigger sprayer to separate it from trigger sprayer
4. Align filter end of the tubing with line drawing shown in Figure 2
5. Using scissors, cut tubing to the length corresponding with Figure 2
6. Re-insert tubing into trigger sprayer
7. Screw adapter onto bottle
8. Guide the tubing of the trigger sprayer through the opening on the adapter and into the bottle
9. Screw trigger sprayer into adapter
10. The spray adjustment pattern can be varied from mist to stream by turning the spray adjustment nozzle as shown in Figure 1.

Figure 1


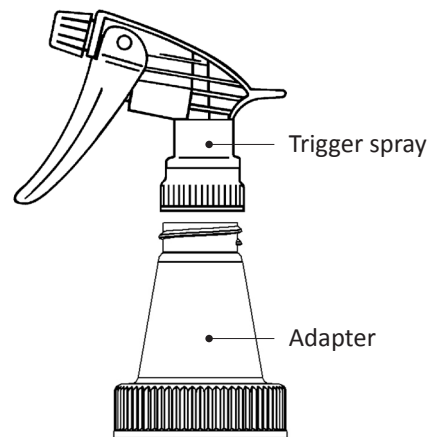


Figure 2



Cut tubing here

Align end of filter with picture

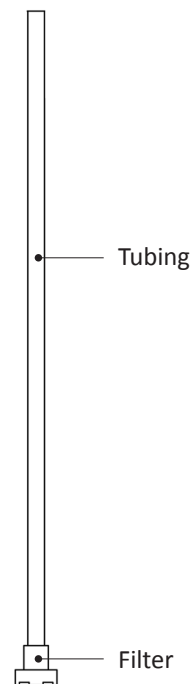

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

