



# Clavies® Biohazard Bag Holders

Catalog No. F13192-0001, F13192-0002, F13192-0003

## 1. Small F13192-0001 (7.3 x 14 x 18.5")

### Parts List

1. Tray (1)
2. Wire Frames (2)
3. Plastic Wire Anchors (4)

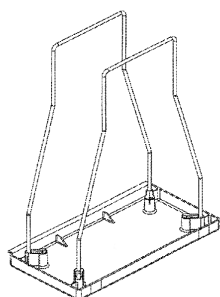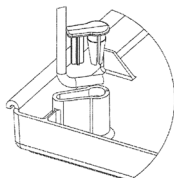

- a. Turn the tray (Part #1) right side up so the four wells are visible.
- b. Push the ends of the stainless steel wire frames (Part #2) into the wells.
- c. Push the four plastic wire anchors (Part #3) into the wells until they are fully seated on the rim of the well and lock the wire frame into place.
- d. Drape the appropriate size bag over the wire frames and roll the edge over the top of the frames.

**Recommended bags for small Clavies® Biohazard Bag Holder are**  
**SP Bel-Art – F13160-0009, F13160-0005,**  
**H13185-1224 and F13164-1923**

### You May Also be Interested in:

#### Clavies® Biohazard Bag Holder

Improve sanitation by adding a cover to the bag holder.


Catalog No. F13192-0102




## 2. Medium F13192-0002 (14 x 14 x 21.5"), Large 13192-0003 (14 x 14 x 27.2")

### Parts List

1. Trays (2)
2. Wire Frames (2)
3. Plastic Wire Anchors (8)
4. Wire Locking Clips (2)

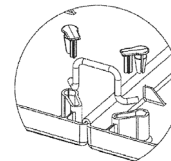

- a. Turn the two trays (Part #1) right side up so the eight wells are visible.
- b. Join the two trays together by slipping the rolled edge of one over the straight edge of the other.
- c. Install the two stainless steel wire locking clips (Part #4) in the center well to lock the trays together.
- d. Push the ends of the stainless steel wire frames (Part #2) into the wells.
- e. Push the eight plastic wire anchors (Part #3) into the eight wells until they are fully seated on the rim of the well and lock the wire frames and wire locking clips into place.
- f. Drape the appropriate size bag over the wire frames and roll the edge over the top of the frames.

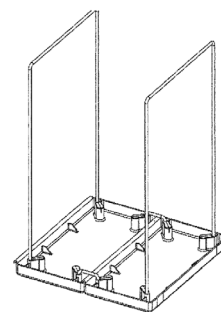

**Recommended bags for medium Clavies® Biohazard Bag Holder are SP Bel-Art F13161-0009, F13161-0005 & H13185-2430**  
**Larger bags can also be substituted.**

**Recommended bags for Large Clavies® Biohazard Bag Holder are SP Bel-Art F13162-0009, F13162-0005, H13185-2436, F13164-2535, F13165-2535, F13164-3138 & F13165-3138**  
**Larger bags can also be substituted.**