


Plastic Economy Support Stand

Catalog No. F18302-0000

To assemble the support stand:

1. Press out the small removable disc located near the edge on the top smooth surface of the base, this will create a mounting hole for the vertical rod.
2. By hand, thread one of the supplied hex nuts onto the threaded end of the stainless steel rod as far as it will go.
3. Slip the threaded end of the rod into the hole in the base and secure it by screwing the second nut onto the threaded rod protruding under the base.
4. Securely tighten the top and underside nuts with a wrench.

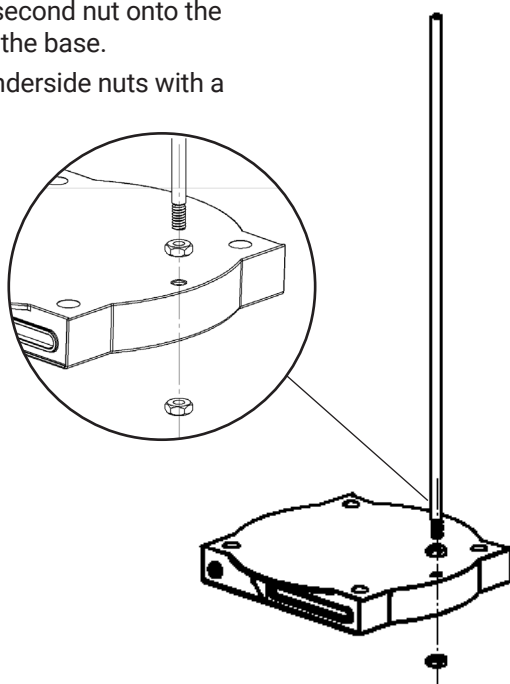

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

