
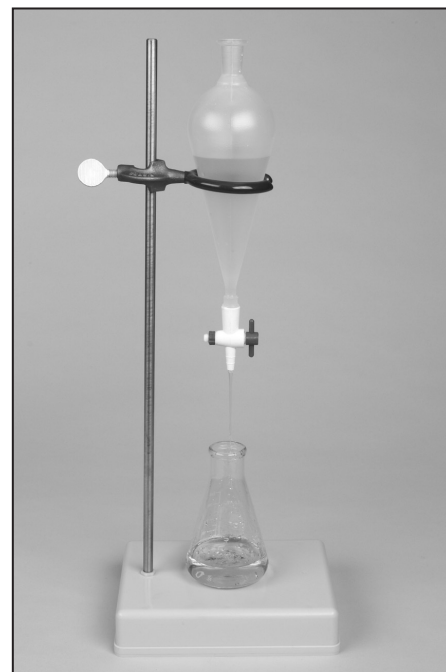


Compact Support Stand

Catalog No. H18303-0005

Directions for use:

1. Screw the hex nut onto the threaded end of the stainless steel rod and continue turning until it stops.
2. Screw the stainless steel rod into the desired hole in the support base and continue turning until it stops.
3. Using pliers, turn the hex nut clockwise until it makes contact with the base. Do not over tighten the hex nut as this may cause it to crack. The purpose of the hex nut is to prevent the rod from turning during use.
4. Screw the hex plug into the unused hole with a 3/16" Allen key until the top of the plug is even with the surface of the polypropylene base.

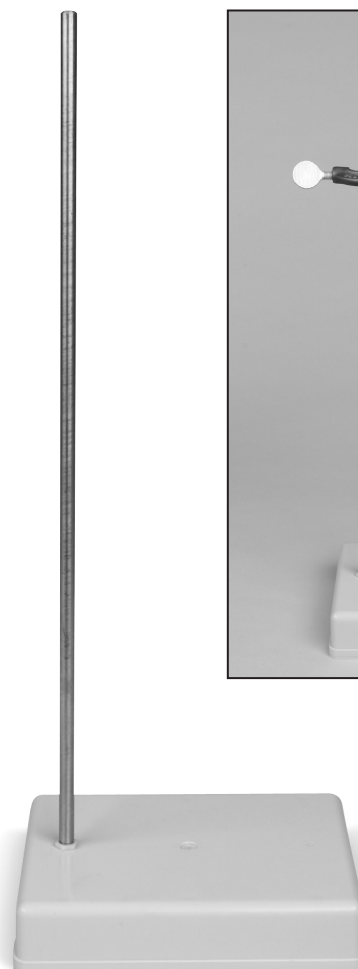

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.



918303038 – Rev 4