

# The Collector™ Blood Tray


Catalog No. F18631-0715

## Locking clips for handle

Make sure both LOCKING CLIPS are installed before using the Blood Tray with the handle

### To install handle LOCKING CLIPS:

- 1) Fit the round posts at the bottom of the handle through the mounting holes on each end of the Blood Tray.
- 2) Slide the LOCKING CLIPS completely down along the inside wall surface of the Blood Tray as shown in FIGURE 1 until the round posts of the handle are trapped and prevented from sliding back out of the mounting holes.
- 3) Test to make sure that the handle is secured to the Blood Tray

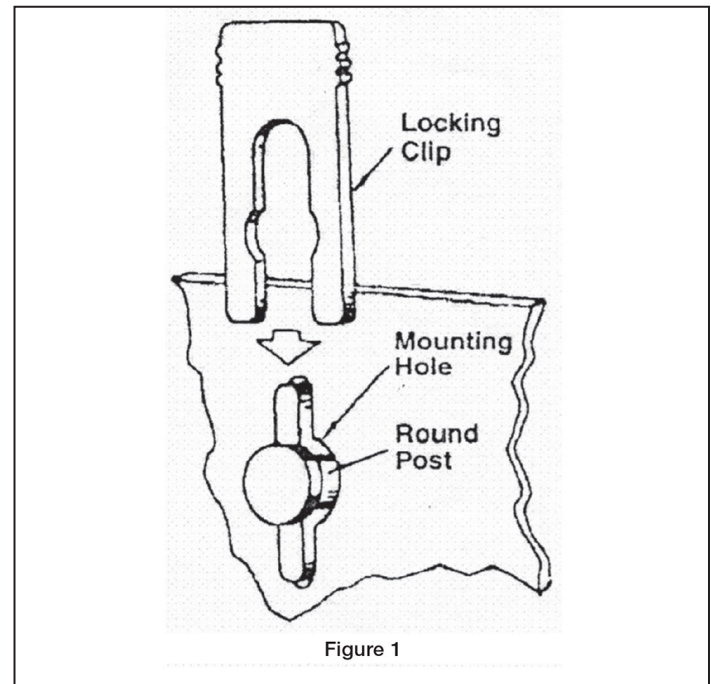

### SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | [cservice@scientificproducts.com](mailto:cservice@scientificproducts.com)

[belart.com](http://belart.com)

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit [scientificproducts.com](http://scientificproducts.com).

