

Switch-Grid Racks

Tube Size	Half Rack Capacity	Full Rack Capacity	Grid Color
10-13mm	42 (7x6)	84 (14x6)	Orange
13-16mm	30 (6x5)	60 (12x5)	Blue
16-20mm	20 (5x4)	40 (10x4)	Green
20-25mm	12 (4x3)	24 (8x3)	Fuchsia
25-30mm	6 (3x2)	12 (6x2)	Yellow


Rack components:

- Switch-Grid Rack Frame
- Grids (4/pack)
- Rubber Feet (4 pieces per rack)

1. Install 4 rubber feet in the corner holes on the bottom of rack by inserting the stems through the bottom of the frame (A). Pull stems until rubber feet pop into place. Excess stem may be trimmed with scissors.
2. Place grid between top and middle shelf of rack with large tabs toward center of rack (B).
3. Insert large tabs onto center rib of frame (C). Compress spring latch (3 small tabs) toward center of rack allowing rack to drop into outer rib (D).
4. Place second grid above top shelf of frame and repeat step 3.
5. Repeat process for opposite side of rack.

Grid Color/Size Chart

Catalog No.	Tube Size	Grid Color
F18745-1300	10-13mm	Orange
F18745-1600	13-16mm	Blue
F18745-2000	16-20mm	Green
F18745-2500	20-25mm	Fuchsia
F18745-3000	25-30mm	Yellow

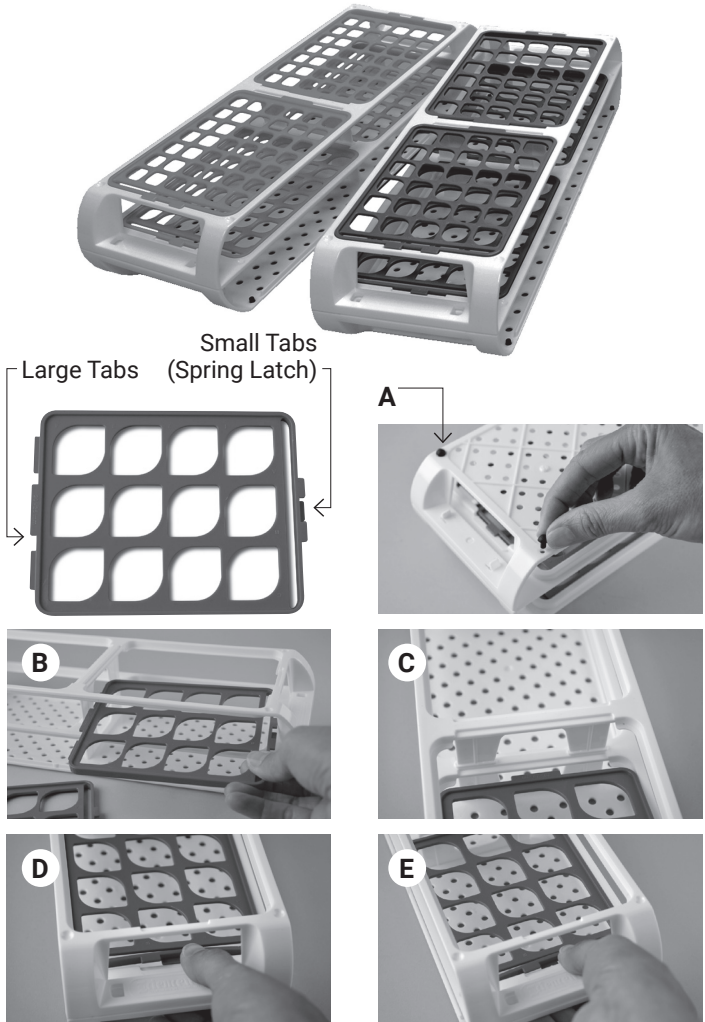

A rack without grids is also available for purchase to make use of your extra Grids.
Look for Catalog Item: F18745-0045



SP Bel-Art
2024 Broad Street, Pocomoke, MD 21851 USA
Tel: 973-694-0500 | cservice@scientificproducts.com
belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

