


No-Wire™ Half Racks

1. Place the base on a flat work surface with the word "BOTTOM" facing down.
2. Bend up one end plate. Snap the bottom tier end tabs into the lower set of slots in the end plate. Be sure the bottom tier imprinted numbers face up.
3. Position the top tier so that imprinted numbers face up and align with the numbers on the bottom tier. Snap the top tier end tabs into the upper set of slots in the end plate.
4. Bend up the opposing end plate and snap the bottom tier and then top tier end tabs into their matching slots in the end plate.
5. Bend up one side tier support and snap the bottom tier and then top tier side tabs into their matching slots in the side tier supports.
6. Bend up the remaining side tier support and snap the bottom tier and then top tier side tabs into their matching slots in the side tier supports.

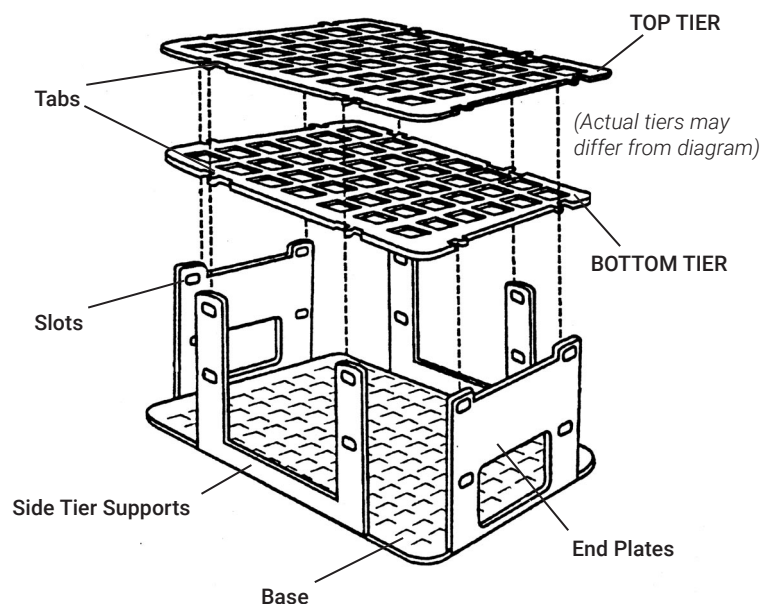

Submersible and works well in water baths

Autoclavable at 250° F (121°C)

Empty racks can be stacked.

Note: Do not use with liquid nitrogen.



SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.



918748001 Rev. 2

No-Wire™ Half Racks

Stackable • Autoclavable • Space Saver • Color Coded • Submersible • End Plates for Labeling • Easy Tube Identification

NO-WIRE™ AUTOCLAVABLE VIAL AND BOTTLE RACKS

Catalog No.	Hole Size mm	Place	Rows	Fits Vials or Bottles
☐ 18514-0013	13.38	42	6 x 7	13 mm
☐ 18514-0016	16.38	30	5 x 6	16 mm
☐ 18514-0020	20.37	20	4 x 5	20 mm
☐ 18514-0025	25.40	12	3 x 4	25 mm
☐ 18514-0030	30.48	6	2 x 3	30 mm

Dimensions: 5.08" x 4.15" x 1.70" high

NO-WIRE™ AUTOCLAVABLE CUVETTE RACKS

Catalog No.	Holds Cuvettes	Place	Rows	—
☐ 18520-0000	10 mm	42	6 x 7	—

Dimensions: 5.08" x 4.15" x 1.70" high

NO-WIRE™ AUTOCLAVABLE MICRO TUBE RACKS

Catalog No.	Holds	Tubes	Place	Rows
☐ 18838-0015	1.5 ml	micro	42	6 x 7
☐ 18838-0005	.5 ml	micro	63	7 x 9
☐ 18838-0400	400µl	micro	80	8 x 10

Dimensions: 5.08" x 4.15" x 1.70" high

NO-WIRE™ AUTOCLAVABLE MICRO TUBE RACKS

Catalog No.	Holds Tubes	Place	Rows	—
☐ 18748-0013	13 mm	42	6 x 7	—
☐ 18748-0016	16 mm	30	5 x 6	—
☐ 18748-0020	20 mm	20	4 x 5	—
☐ 18748-0025	25 mm	12	3 x 4	—
☐ 18748-0030	30 mm	6	2 x 3	—

Dimensions: 5.08" x 4.15" x 2.83" high

NO-WIRE™

