


No-Wire™ Grip Racks

Stackable • Autoclavable • Space Saver • Retains Tubes for Decanting • Submersible • End Plates for Labeling • Easy Tube Identification


No-Wire™ Grip Test Tube Racks are perfect for every laboratory use, with distinct advantages over other plastic racks. All test tubes are held securely in place, which is especially useful in decanting the content of all tubes simultaneously. Now, with the addition of this economical new line of No-Wire™ Grip Racks, SP Bel-Art provides you with the widest variety of test tube racks.

Catalog No.	Holds Tubes	Places	Rows
18749-0000	10-13mm	90	6 x 15
18749-0001	15-16mm	60	5 x 12
18749-0002	18-20mm	40	4 x 10

Assembly instructions on reverse side.

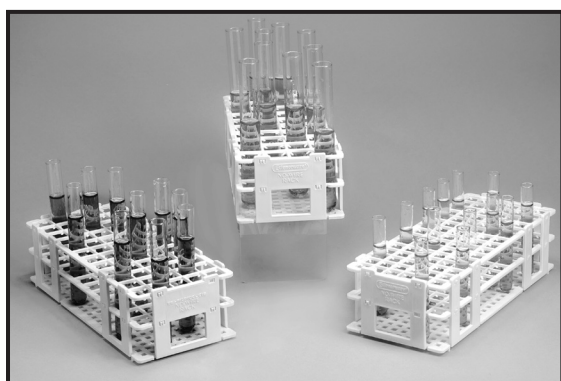

No-Wire™ Grip Racks

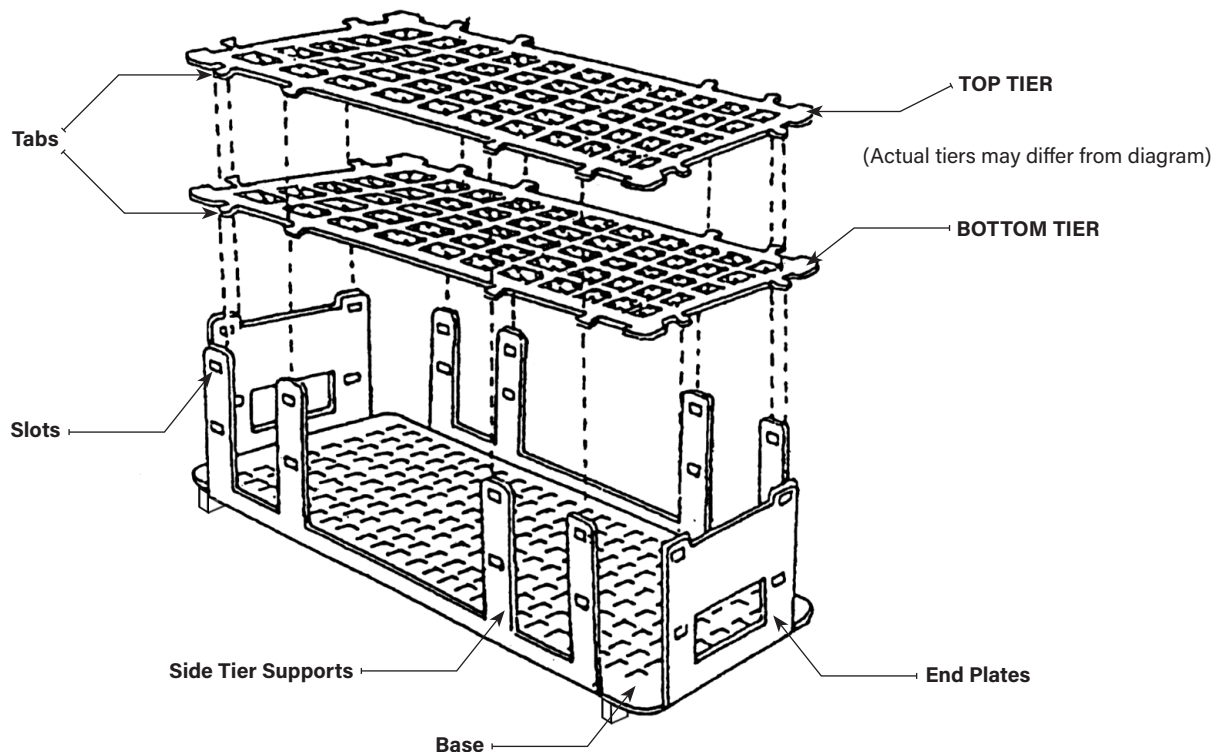

1. Place the base on a flat work surface with the word "BOTTOM" facing down.
2. Bend up one end plate. Snap the bottom tier end tabs into the lower set of slots in the end plate. Be sure the bottom tier imprinted numbers face up.
3. Position the top tier so that imprinted numbers face up and align with the numbers on the bottom tier. Snap the top tier end tabs into the upper set of slots in the end plate.
4. Bend up the opposing end plate and snap the bottom tier and then top tier end tabs into their matching slots in the end plate.
5. Bend up one side tier support and snap the bottom tier and then top tier side tabs into their matching slots in the side tier supports.
6. Bend up the remaining side tier support and snap the bottom tier and then top tier side tabs into their matching slots in the side tier supports.

Submersible and works well in water baths

Autoclavable at 250° F (121°C)

Empty racks can be stacked.

Note: Do not use with liquid nitrogen.



SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

