



Test Tube Capper

Catalog No. F18799-0000 (Patent Pending)

- Eliminate fatigue and soreness of hands and fingers from repetitive task of applying test tube caps
- Fits in the palm of your hand
- Designed for use with recessed center test tube caps

Instructions for use

1. Grasp capper with the black top resting in the palm of your hand.
2. Place cap on the test tube.
3. Using the Test Tube Capper, gently press cap securely onto the test tube.

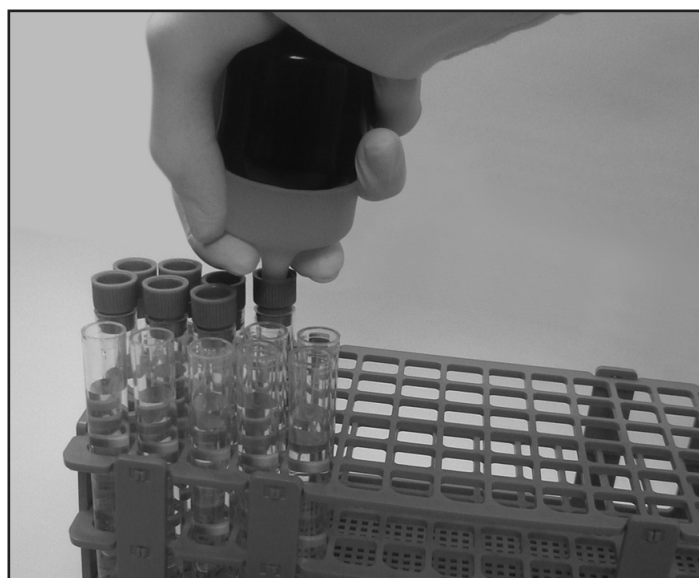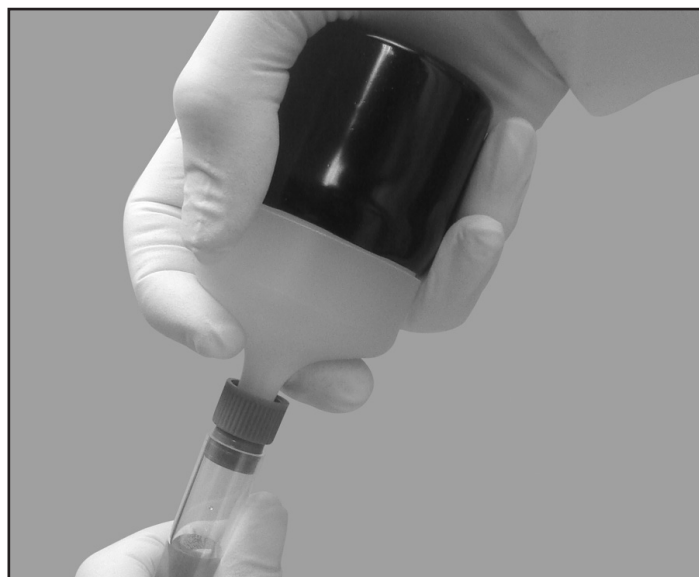

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.



918799002 Rev 2