



Rack Slanter

Catalog No. F18890-0100


1. Locate the word "BACK" on the board and position the board so that the word "BACK" is facing the benchtop (at selected angle).
2. Locate the board supports that are labeled with the letter "L" for the left support and the letter "R" for the right support.
3. Stand the left and right rack support on the benchtop with enough distance in between each to fit the backboard into the desired support slots.
4. Insert the end of the backboard containing a slot "L/B" (**Left Board**), into one of the desired support slots labeled "L/S" (**Left Support**) corresponding to the desired angle.
5. Insert the end of the backboard containing a slot "R/B" (**Right Board**), into the desired support slots labeled "R/S" (**Right Support**) corresponding to the desired angle.
6. Once the board is at the desired angle, insert the four peg supports as shown in the diagram.
7. Place the laboratory test tube rack on the peg supports.

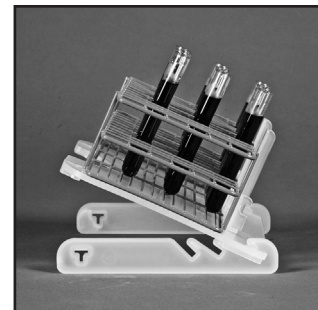


25°

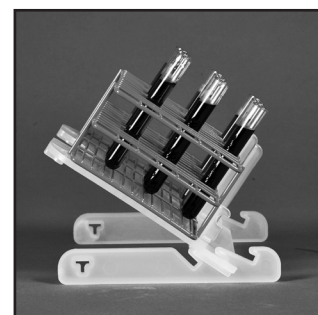

35°

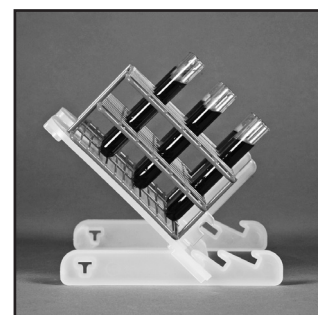

45°

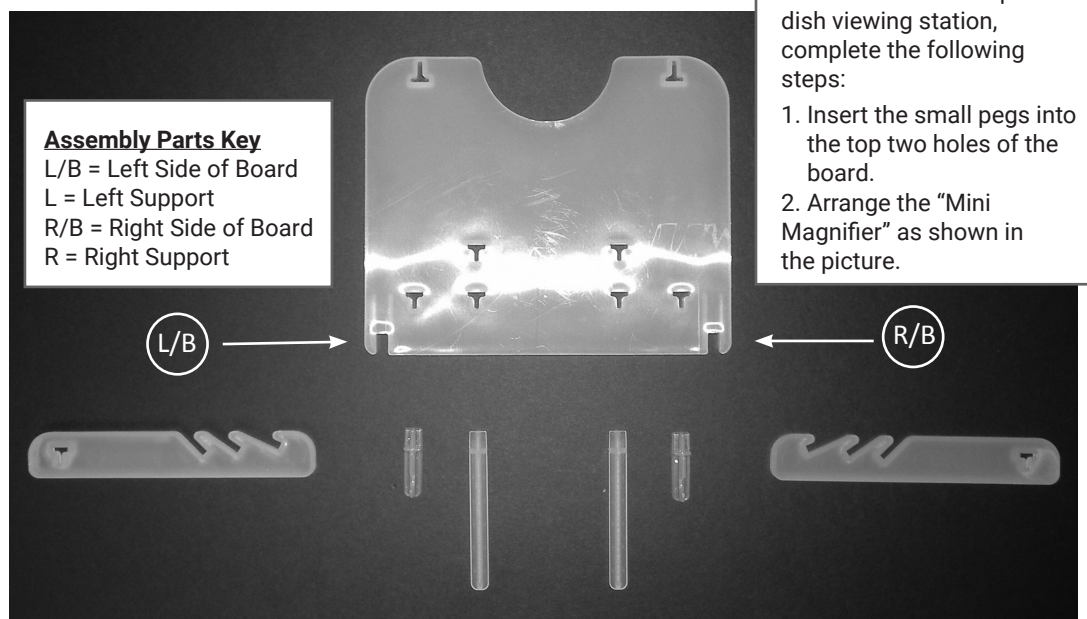

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.



918890106 – Rev 3