

Lab-Aire® II Wall Mounted Unit Labware Drying Rack

Catalog No. 18933-0011 • 18933-0013 • 18933-0014


Wall Mounting:

2 hooks or 2 #8 screws should be horizontally fastened on to a vertical mounting surface exactly 10 inches apart. In the case when screws are used, make sure the screw heads are backed off 3/16 inch from the mounting surface to allow for the thickness of the key slots of the hanging tabs.


Peg Assemblies:

On each panel, pegs should be arranged with the small pegs on the top rows and large pegs on the bottom rows. To ensure pegs are not impacted too tightly for removal or rearrangement, do not force the pegs into the peg holes.

Description	Catalog No.
Bag of 19 Small Pegs & 6 Plugs	18931-0001
Bag of 19 Medium Pegs & 6 Plugs	18931-0002
Bag of 19 Large Pegs & 6 Plugs	18931-0003

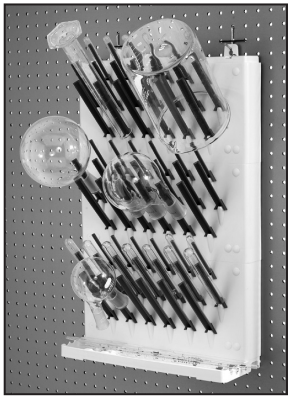

18933-0013

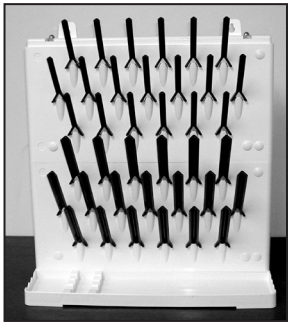

18933-0014



SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA
Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

