


PiRack™ Pipettor Holder

Catalog No. 18954-0000 Triple (clamp and 2 extenders) A&B
18954-0002 Extenders (4 extenders only) B only

The SP Bel-Art PiRack™ Pipettor Holder is designed to hold many major brands of pipettors, facilitating easy to reach shelf or bench top placement and creating flexibility in locating and storing your pipettors. The SP Bel-Art PiRack™ clamps and extenders can be assembled to hold virtually any number of pipettors. The only limitation is the space in your lab.

Note: It is recommended that for every clamp employed, a maximum of two extenders be used to obtain optimum load distribution.

Assembly:

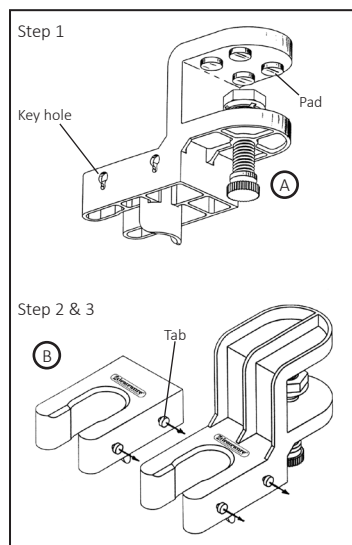


1. Hold all components with the brand logo facing up.
2. Apply the four (4) non-skid adhesive backed pads evenly to the underside of the clamp's upper lip.
3. Align the round assembly tabs on the side of the extender with the corresponding key holes in the side of the clamp.
4. Insert the tabs into the holes. Hold upright. Position your thumbs on the extender and apply an even, downward force

until the tabs lock into place. If necessary, adjust the extender until its top surface is even with the surface of the clamp.

5. To assemble the other extender to the clamp, repeat instructions 3 through 4.

Note: The tabs will now be on the side of the clamp and the holes will be in the side of the extender.

Shelf Attachment:

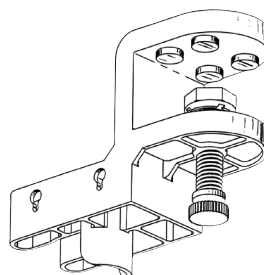

Note: The hex nut can be removed from the thumbscrew for wider shelves.

1. Position the shelf edge between the upper lip of the clamp and hex nut on end of locking thumbscrew.
2. Position clamp until the edge of the shelf is located squarely in the "U" shape of the clamp.

Note: A minimum shelf width of 2 3/16" is required.

3. Tighten thumbscrew until PiRack™ Pipettor Holder is securely fastened.

Note: The maximum shelf thickness is 1 1/2".



SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

