


Drying Tubes With Tube Fittings

Catalog No. 19960-0000, 19961-0000, 19962-0000

Extra Fittings Cat. No. 19965-0000

The drying tubes are made of polyethylene and should not be placed in an autoclave. The end fittings are molded of polypropylene and are autoclavable at 250°F (121°C). The tubes are supplied empty and may be filled with a drying agent. The drying tubes may vary slightly in diameter. For this reason, the end fittings are molded with a “stepped” inner diameter. If the tubes are on the upper limit of the diameter ranges, the fitting will slip on only about 1/4” (6.4mm). It is not necessary for the tube to enter the fitting completely. If the tubes are on the lower limit of the diameter range, they will fully enter the fitting.

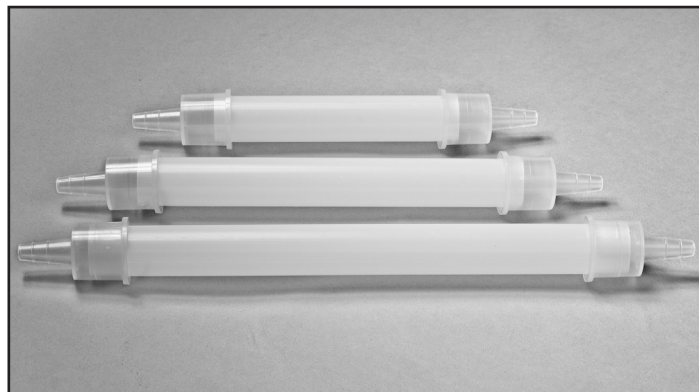

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

