

Splash Shields & Weighted Base Bracket

Catalog No. F24966-0008, F24966-0006, F24966-0005 & F24966-0007

Splash Shields


SINGLE ARM SHIELD

1 arm with neck fitting
2 handles with bolts/washers
(these shields need
1 attachment bracket!)

TWIN ARM SHIELD

2 arms with neck fittings
2 handles with bolts/washers
(these shields need
2 attachment brackets!)

PARTS LIST/PLEXIGLAS® ASSEMBLY

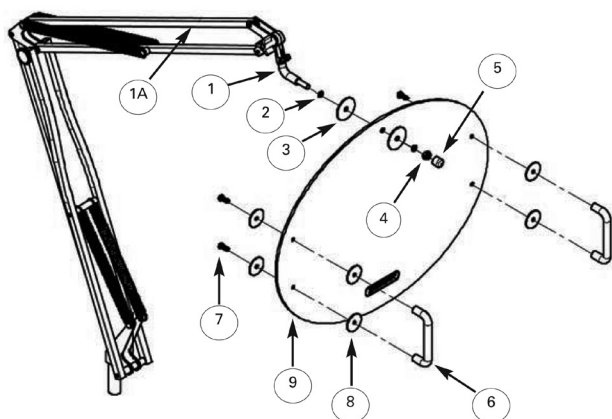

PART NUMBER AND QUANTITY

- | | |
|---------------------------|-------------------------|
| 1A. Arm (1). | 5. Cap, Red (1). |
| 1. Neck (1). | 6. Handle (2). |
| 2. Washer, Lock (2). | 7. Screw (4). |
| 3. Washer, 1.25 O.D. (2). | 8. Washer, 1" O.D. (8). |
| 4. Lock Nut (1). | 9. Shield, Oval (1). |

NOTE: ALL SHIELDS NEED ATTACHMENT BRACKETS!

Weighted Base
Adapter



F24966-0008

Fixed Table
Bracket



F24966-0006

C
Clamp




F24966-0005

Wall
Bracket



F24966-0007

WEIGHTED BASE BRACKET

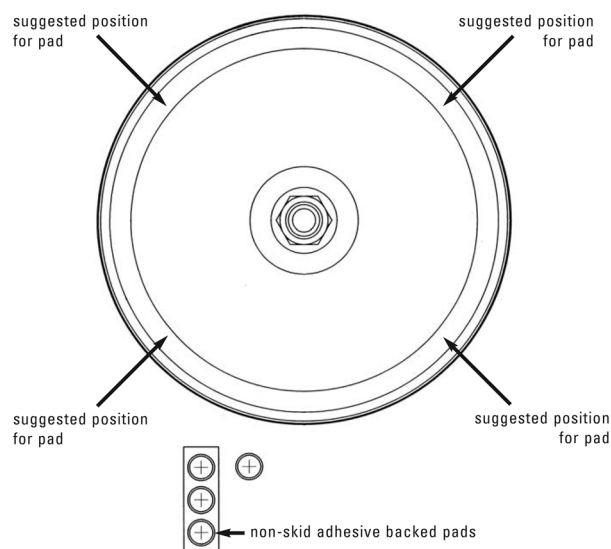

1. Apply the four (4) non-skid adhesive backed pads, included with your SP Bel-Art Weighted Base Bracket, to the four positions on the underside of the base.
2. Insert the shield arm into the receptor on the upper side of the weighted base.



SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.



924966107 Rev. 2