

## Buchner Funnel Trowel

Catalog No. F36820-0118


The SP Bel-Art Buchner Funnel Trowel is designed to smooth the cake (i.e. the precipitate) in 18, 24 and 36" Buchner Funnels, allowing for more efficient filtration.

### Directions for Use:

Begin filtration in the Buchner funnel. Holding the SP Bel-Art Buchner Funnel Trowel at a slight angle in either left or right hand, sweep the trowel over the surface of the cake and smooth over any noticeable cracks or fissures. Continue smoothing the cake as needed until desired filtration level is reached.

### Care and Cleaning:

- Do NOT autoclave product.
- Compatible with industrial dishwasher water temperature and cleaners
- Recommended detergent: Aquet®, Catalog no. F17094-0030, 1 gallon

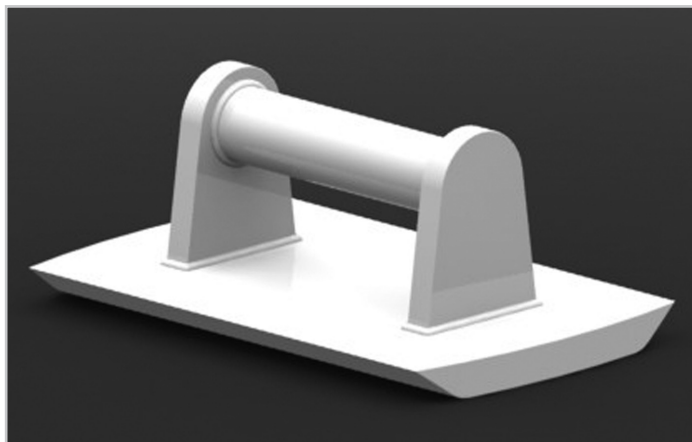

### SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | [cservice@scientificproducts.com](mailto:cservice@scientificproducts.com)

[belart.com](http://belart.com)

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit [scientificproducts.com](http://scientificproducts.com).

