
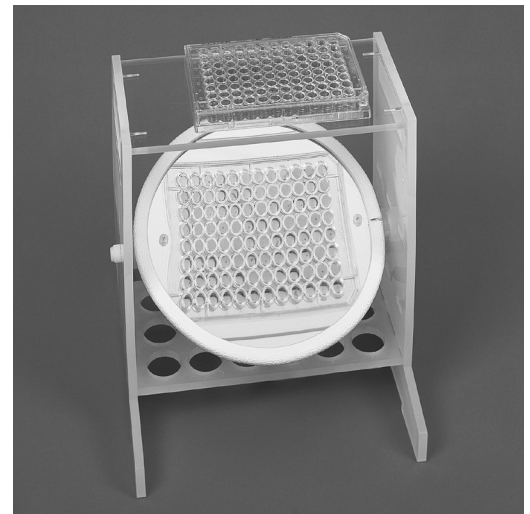


Plate Reader

Catalog No. 37861-0000

1. Place the reader on a flat, level surface.
2. Place a petri dish with cell colony growth on the top clear acrylic plate. 96 well microliter plates can also be used.
3. Loosen the two slide nuts on the frame to adjust the mirror. Adjust the mirror horizontally and tilt the mirror to reflect the cells under observation. Tighten the nuts to hold the mirror in place.
4. Both the mirror surface and the top plate should be kept clean for best observation results. Do not use abrasive cleaners or any type of scrub pad which might scratch the acrylic surfaces. Cleanware® Lab Equipment Wipes (Bel-Art catalog No. F17086-0000) leave a clean, non-smear surface. A dampened, clean, soft cloth with mild detergent may also be used. Remove any residual detergent with the damp cloth. Dry with a clean, soft cloth.
5. Store the reader where it will not be subject to dust. A plastic cover or bag will keep dust out.

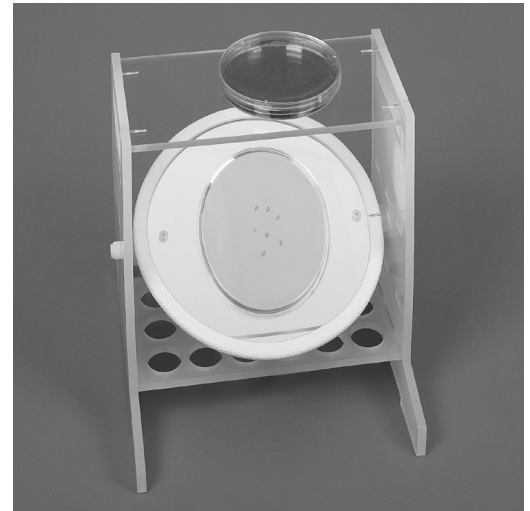

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

