


I N S T R U C T I O N S



Bel-Blotter™ Catalog No. 37876-0002

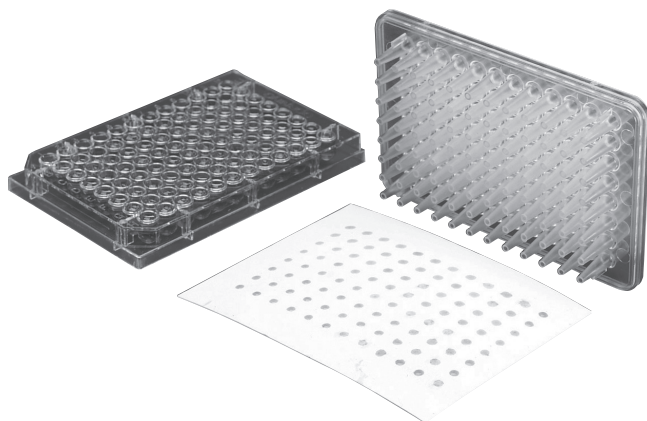

1. The Bel-Blotter, as supplied, is non-sterile, however, it is molded of autoclavable polycarbonate and should be sterilized before use.
2. This replicating tool will do the tedious replicating tasks with little effort and increased productivity.
3. The Bel-Blotter will fit all types of 96 well plates ranging from flat, vee or round bottom plates to the 0.2ml thin walled PCR plates and tubes.
4. By simply placing the Bel-Blotter into the plate, the tips will pick up fluid, up to 10µl per tip.
5. The filled Bel-Blotter is then placed into the receiving medium, whether it be the wells of plates, hybridization membranes or agar.
6. The resulting products of the Bel-Blotter can be used for replication of recombinant DNA libraries, inoculation of filters for colony hybridization, phage typing and several other applications.

Let others know what you think. Write a product review at www.belart.com

661 Route 23 South, Wayne, NJ 07470

TEL: 973-694-0500 • FAX: 973-694-7199 • www.belart.com

Bel-Art – SP Scienceware assumes no obligation or liability for any advice furnished by it, or for results obtained with respect to these products.

All such advice is given and accepted at the buyer's risk.

© 2020. All Rights Reserved.

Bel-Art® – SP Scienceware® is a division of SP Industries, Inc.



937876201 – Rev. 2



YOU MAY ALSO LIKE

Bel-Art Hand-Held Colony Counter

A cost effective and accurate way to count colonies of bacteria or other microorganisms grown on an agar plate. This simple to operate device marks, counts, and gives an audible and visual alert each time a count is registered. The lightweight, pen-style design reduces the likelihood of error, combining all aspects of counting into one motion. F37862-0002 is CE marked.

- Fast, accurate and economical sample counting and sample count verification
- Stores up to 35 sample counts
- Verify sample counts with up/down counting
- Displays sum total of all colonies counted
- Backlit display for easy viewing in low light conditions
- Audible alert can be turned off if desired
- Counting function can be locked to prevent error during breaks between counting
- Long life, red and black felt-tip markers (included) work on plastic or glass and are oil and water resistant; ink is removable with ethyl alcohol
- Includes 3V lithium battery and protective storage case



Catalog No.	Description	Quantity / Pack
F37862-0002	Hand-Held Colony Counter, CE marked	1

For more information, visit www.belart.com