

ProCulture® Round Microcentrifuge Floating Bubble Rack

Catalog No. F18875-1600, F18875-1000, F18875-0400

Use as follows:

To Fill:

1. Place the rack on a bench top.
2. Insert the empty tubes.
3. Fill tubes.
4. Close caps securely
5. Lift rack slightly & push down on tubes to lock securely in rack.
6. The rack is now ready for flotation

To Empty:

1. Place filled rack on the bench.
2. Push down the rack to release tubes.
3. Remove tubes.




You May Also be Interested in: Microcentrifuge Tube Refrigerator Rack *Keep Your Microtubes Cool!*

- Two Piece Styrene Tray/Rack
- Bottom Section Designed to Hold Water or Ice for Stabilizing Temperatures
- Large Capacity Stackable Racks
- Available in Two Sizes to Hold a variety of 400µl or 1.5ml tubes.

Catalog No. H18901-0000

H18907-0000

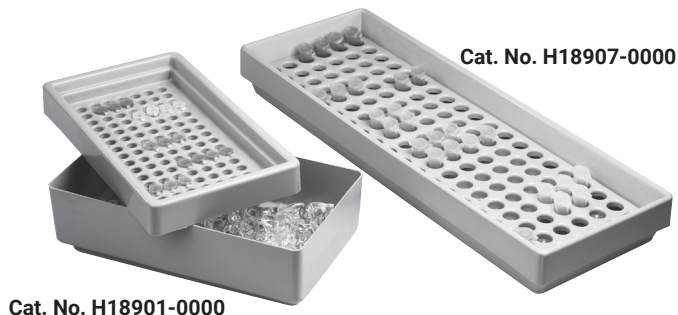

Cat. No. H18901-0000

Cat. No. H18907-0000

Your opinion counts! Write a product review at www.scientificproducts.com