


ProCulture® 60mm Petri Dish Rack

Catalog No. F18991-0060

Assembly:

- 1. Refer to drawing for orientation of parts.
- 2. Assemble the center divider to the base securing with one nut on each of the two threaded posts.
- 3. Assemble each of the rods to the base, secure each rod with one nut on the threaded post.

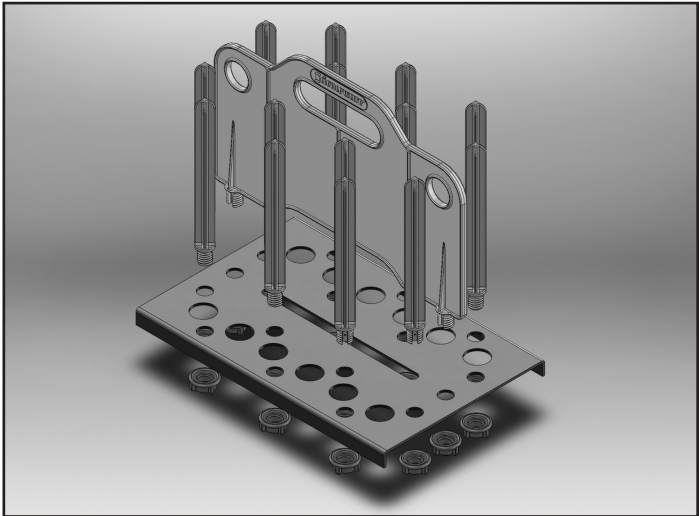


You May Also be Interested in:

Write-On™ Tape Dispensers make labeling easy and fast


Producing legible labels is so much easier on a smooth, flat platform than on curved surfaces such as bottles or vials. The new SP Bel-Art Write-On™ Tape Dispenser provides this convenience by always having labeling tape positioned on the large writing platform, ready to make into a label. Select either the benchtop model with weighted base or economy model that can be wall mounted to save space. All models are easy to use by both left and right handed users.

Each dispenser can hold 5 or more tape rolls with widths from 12 - 25mm (1/2 - 1"). Ideal for use with the full line of SP Bel-Art, Write-On™ Tapes, F13463 series, available in 6 colors and three widths, 12, 19, and 25mm (1/2", 3/4" and 1"). This advanced tape adheres to any clean surface and resists moisture and most solvents. It withstands long exposure to temperatures as low as -73°C (-68°F) and up to 80°C (176°F) yet peels off without leaving a sticky residue.

Contact your preferred distributor to order these and other SP Bel-Art products.

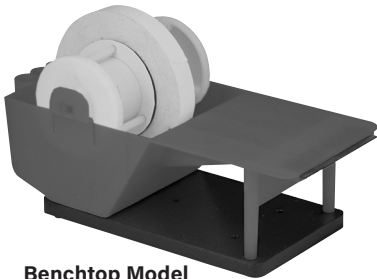

Benchtop Model
Includes weighted base

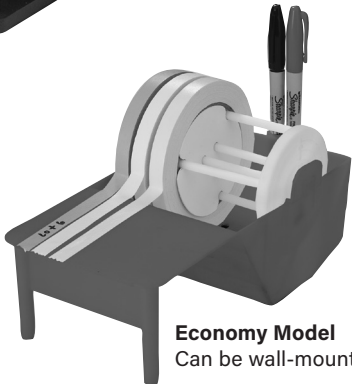

Economy Model
Can be wall-mounted

Catalog No.	Description	Dimensions
F13472-0000	Benchtop Model	302 x 143 x 121 mm H (12 x 5 ⁵ / ₈ x 4 ³ / ₄ " H)
F13473-0000	Economy Model	302 x 143 x 106 mm H (12 x 5 ⁵ / ₈ x 4 ³ / ₁₆ " H)



SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA
Tel: 973-694-0500 | cservice@scientificproducts.com
belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.

