

Aspirator Pump

Catalog No. F32947-0000

The Bel-Art® Aspirator Pump is molded of high-density polyethylene, is equipped with an integral check valve, and removable tube connector. The tube connector may remain attached to the flexible tubing leading from your filtering flask.

To reconnect the flask, simply plug the connector into the aspirator. An enclosed length of 3/8" O.D. polyethylene tubing is provided to insert into the discharge end of vacuum pump to reduce splash.

The Bel-Art® Aspirator Pump operates efficiently on water pressures of 11 psi (75842 N/m²) and above. Maximum attainable vacuum is 27.4" Hg (696mm) with an approximate water flow of 1.5 gallons per minute (6 liters).

The Aspirator Pump is excellent for use with corrosive filtrates and may also be used as a vacuum suction tool to draw up spilled acids.



You May Also be Interested in:

"SPACE SAVER" VACUUM DESICCATORS

Choice of Three Sizes and Two Models

A. Polypropylene Bottom Desiccators



42010-0000

42020-0000

42025-0000

B. Polycarbonate Bottom Desiccators



42012-0000

42022-0000

42027-0000

LARGE INTERIOR VOLUME • EXTRA HEAVY WALLS


Catalog No.	Flange O.D.	Inside Diameter	Plate	Overall Height	Clearance Above Plate
A. Polypropylene Bottom					
42010-0000	171mm (6 ³ / ₄ "	149mm (5 ⁷ / ₈ "	140mm No.3	206mm (8 ¹ / ₈ "	121mm (4 ³ / ₄ "
42020-0000	230mm (9 ¹ / ₈ "	197mm (7 ³ / ₄ "	190mm No.4	260mm (10 ¹ / ₄ "	157mm (6 ¹ / ₈ "
42025-0000	273mm (10 ³ / ₄ "	240mm (9 ⁵ / ₈ "	230mm No.5	311mm (12 ¹ / ₈ "	197mm (7 ³ / ₄ "
B. Clear Polycarbonate Bottom					
42012-0000	171mm (6 ³ / ₄ "	149mm (5 ⁷ / ₈ "	140mm No.3	206mm (8 ¹ / ₈ "	121mm (4 ³ / ₄ "
42022-0000	230mm (9 ¹ / ₈ "	197mm (7 ³ / ₄ "	190mm No.4	260mm (10 ¹ / ₄ "	157mm (6 ¹ / ₈ "
42027-0000	273mm (10 ³ / ₄ "	240mm (9 ⁵ / ₈ "	230mm No.5	311mm (12 ¹ / ₈ "	197mm (7 ³ / ₄ "

Vacuum-Pressure Indicator

Catalog No. 42042-0000

Attaches Easily

The Vacuum-Pressure indicator fits between the desiccator dome and the desiccator stopcock.

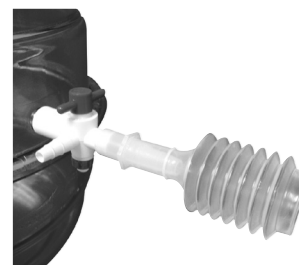

SP Bel-Art

2024 Broad Street, Pocomoke, MD 21851 USA

Tel: 973-694-0500 | cservice@scientificproducts.com

belart.com

Individual use varies. SP Bel-Art assumes no obligation or liability for any advice furnished by it or for results obtained with respect to these products. All advice is given and accepted at the user's risk. For warranty information and to register your product visit scientificproducts.com.



932947003 – Rev 4